



# EVA WELL

A new vision of women's health

EVA WELL is an essential tool for gynecologic examinations, enhancing visualization and documenting efficacy and ensuring safety.

# EVA WELL

## The essential tool for gynecologic examinations

EVA WELL is a user-friendly, portable, digital visualization tool for vaginal, vulvar and cervical examinations. EVA WELL couples advanced hardware with integrated software to enhance visualization and monitor the safety and efficacy of gynecologic procedures.



### Visualize and capture

Achieve high-quality medical-grade images and video with EVA's 16x optical magnification.



### Collaborate remotely

Consult with peers or remote experts via secure chat and live video case-sharing in the EVA System Online portal.




### Streamline safety and efficacy

Document any abnormalities, add annotations and locations to images, and track patient progress over time.



### Educate and reassure your patients

Share images with your patients. Improve patient communication and information about their procedure and their health.

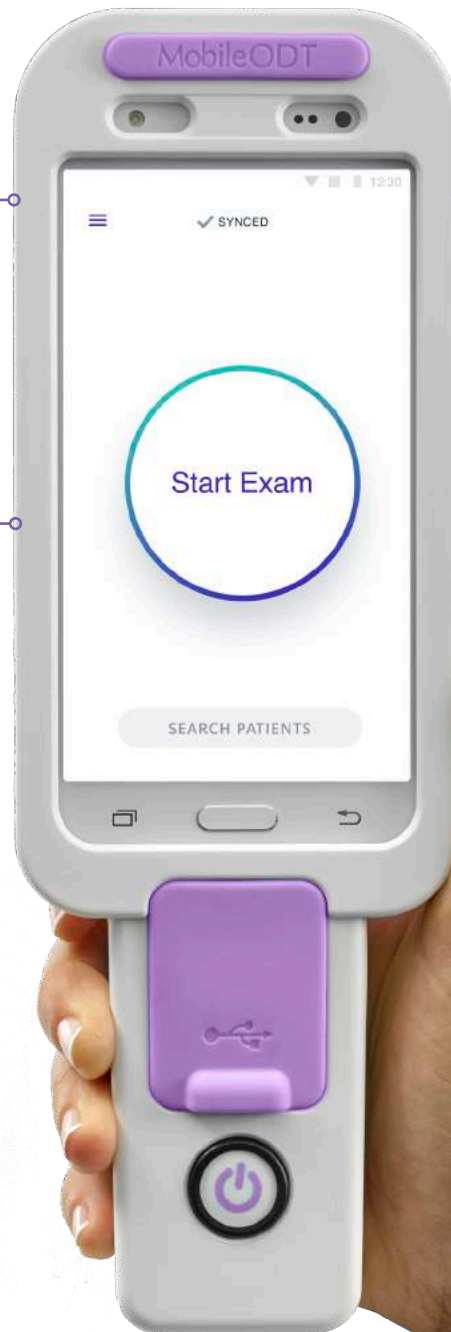

*"With EVA WELL, every vaginal examination I perform is documented. My notes are added to the images and recorded on a secure online portal. I can see the progression over time. I can show my patients their own anatomy and even email them the images after the exam."*

**Dr. Pablo Naranjo**  
Medical Director of the Elite Laser  
Clinic, Madrid



## State of the art medical-grade optics

The EVA System is designed to provide clinicians with high quality images to aid visualization and diagnosis:



### 16 times magnifying lens

EVA WELL has a combined 16x image magnification (4x optical and 4x digital). Its 13 MP camera resolution and 250-400mm working distance allow for a superior image and video capture.

### Crossed-polarizers for glare reduction

EVA WELL implements a reflective wire-grid polarizer to reduce glare. The polarizers allow diffusely reflected light to reach the camera but reject specular glare reflected off the tissue surface, making the images clean and clear.

### LED light

EVA WELL incorporates a white LED light with a 6000k color temperature and polarizer to ensure consistent well-lit visualization and image capture.

### Green filter

EVA WELL's capture includes a digital green filter that has been designed to optimize the contrast of vascular structures and blood.

### Contactless image capture

With the EVA System's hand-waving feature, clinicians are able to capture images with a movement of their hand in front of the display without needing to touch the device itself, maximizing image quality while reducing the likelihood of contamination.

*"Any physician who needs to visualize the vagina, vulva or cervix should be doing it with EVA WELL. Even for exams that could be performed with the naked eye, EVA WELL lets you document your observations and show patients exactly what is happening."*

**Dr. Jonah Mink**

**Medical Director, MobileODT, Philadelphia**

# A solution for every need

At MobileODT, we recognize that a medical platform needs to be as flexible as its customers are unique. And that's a challenge we readily embrace. So whether you're a women's health clinician working independently, or a hospital looking to optimize ROI and services, MobileODT has the solution to get you on your way.

## TECHNICAL SPECS EVA 3.0 PLUS

<b>Operating System</b>	Android 7
<b>Mobile Device</b>	Samsung Galaxy J5 PRO
<b>Processor</b>	1.6GHz Octa-core
<b>Display</b>	Super AMOLED 720x1280 pixels HD 5.2' (131.8 mm)
<b>Camera</b>	- CMOS sensor, 13 MP resolution - 250mm - 400mm working distance - 4.0x optical / 16x with digital zoom
<b>Video</b>	- FHD (1920x1080@30fps) - Playing format: MP4
<b>Connectivity</b>	- 2G GSM, 3G WCDMA, 4G LTE FDD, 4G LTE TDD - Wi-Fi: 802.11 a/b/g/n/ac 2.4+5GHz - Bluetooth V4.1
<b>Power</b>	- Micro USB charger 5V DC - Two batteries: 3000 mAh (mobile device) 6500 mAh (optical unit)
<b>Illumination</b>	- 6000k 3W LED (3.6V) - Wire-grid polarizer for glare reduction
<b>Dimensions and weight</b>	- Colposcope (HxWxD) 24.39 cm x 8.1 cm x 10.41 cm - Packaging (HXWXD) 16cm x 29.4cm x 24.4cm - Weight: 605 gr



### To learn more about EVA WELL contact:

[hello@mobileodt.com](mailto:hello@mobileodt.com)

+ 1 929 376 0061

[www.mobileodt.com](http://www.mobileodt.com)

© 2019 by MobileODT Ltd. All rights reserved. All MobileODT Ltd. products and services mentioned herein are registered trademarks of MobileODT Ltd. All other product, services or research mentioned belong to or are the trademarks of their respective owners. These materials are subject to change without notice. These materials and the data contained are provided by MobileODT Ltd. and its partners for informational purposes only, without representation or warranty of any kind, and MobileODT Ltd. shall not be liable for errors or omissions in this document, which is meant for public promotional purposes. MobileODT reserves the right to deny approval of external usage and reproduction for any reason.